

# CoolTex™ Ag Hybrid Anterior Closure Knee with ROM Hinge

Hybrid design provides the compression of a sleeve below the knee and wrap-style closure on the top of the thigh for improved fit



**CORFLEX®**

Designed for life.


CORFLEXGLOBAL.COM



16"  
Length  
Shown

The CoolTex Ag hinged knee brace offers comfortable support with the added convenience of easy application. Brace features bilateral range of motion hinges delivering full flexion and extension control with settings at nine ROM options. Top finger loops allow for simple pull-on application and adjustment during wear while the release tabs on the hook closure aid in brace removal.

Unique hybrid design features an extended wrap-style closure on the top of the thigh to allow greater adjustability in fit while the lower calf section utilizes a closed sleeve for improved compression. Open patella design offers stability and includes our adjustable Cor-Trak buttress for added tracking control. The brace is manufactured from CoolTex Ag fabric to allow for easy flow of air, transfer of moisture and antimicrobial protection. Lightweight and breathable material allows for greater comfort during extended wear. Available in 13" and 16" lengths in eight sizes from XS through 4XL.

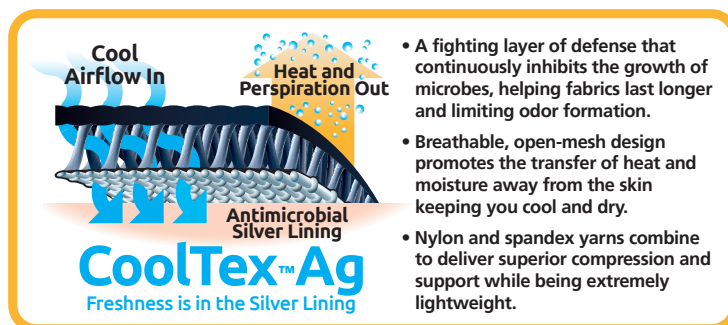

## ADDITIONAL FEATURES:

- Features ROM hinges that allow for full flexion and extension control and a unique design that offers the convenience of a sleeve below the knee and benefits of an adjustable wrap above the knee; available in two lengths
- Indicated for medial/lateral instability, muscle/ligament weakness or injury, ACL/PCL moderate instability, strains or sprains, meniscal injuries and degenerative joint disease
- Manufactured from lightweight, breathable, antimicrobial CoolTex Ag material with open popliteal design; bilateral, aluminum, polycentric hinges with flexion stops at 45°, 60°, 75° and 90°, and extension stops at 0°, 10°, 20°, 30° and 40°; removable, EVA foam buttress



Polycentric Hinge with ROM  
Extension/Flexion Stops

# CoolTex™ Ag 13" and 16" Hybrid Anterior Closure Knee with ROM Hinge



13"



Pull loops assist with application and adjustment during wear.



16"

Open popliteal design and elastic thigh strap promote greater comfort, especially during activities that involve frequent knee flexion.

Thigh closure with extended loop panel allows maximum adjustability during fit.

Elastic strap sections allow for greater comfort and changes in muscle thigh and shape.

Product Number	XS	S	M	L	XL	2XL	3XL	4XL
13"	85-8322	85-8323	85-8324	85-8325	85-8326	85-8327	85-8328	85-8329
16"	85-8342	85-8343	85-8344	85-8345	85-8346	85-8347	85-8348	85-8349

### Sizing Guide

6" Above Mid-Patella	13"-16"	14"-18"	15"-20"	16"-23"	17"-25"	18"-26½"	19"-28½"	21"-30½"
Mid-Patella	11"-12"	12"-13"	13"-14"	14"-16"	16"-18"	18"-20"	20"-22"	22"-24"
6" Below Mid-Patella	11"-12"	12"-13"	13"-14"	14"-15"	15"-16"	16"-17"	17"-18"	18"-19"

### Accessories

Cor-Trak Buttress 10/Pk	88-9000
Cor-Trak Buttress	88-9001

### Suggested Billing Codes

L-1832 (Custom-Fit) or L-1833 (OTS) KO, adjustable knee joints



CORFLEX GLOBAL | 669 EAST INDUSTRIAL PARK DRIVE | MANCHESTER, NH 03109  
CORFLEXGLOBAL.COM | TOLL FREE: 800 426 7353 | FAX: 603 623 4111