

Solar OA™ Osteoarthritis Knee Brace



CORFLEX®

Designed for life.

CORFLEXGLOBAL.COM

Single-visit, custom molding process delivers superior fit with a universal design.




The Solar OA brace is a state-of-the-art solution to the growing demand for custom orthotics. Featuring the SolarTech molding & curing process for near-native zero offloading, a 90-second UV light bath solidifies the brace and ensures custom definition and durability. With a hardness value of D75, once cured and custom-fit, the Solar OA is comparable on the hardness scale to a construction helmet.


Built with a stainless steel hinge, Solar OA leads the market in weight, durability and stabilization. A low-profile hinge provides a superior fit for the patient while silicone cuffs bonded to malleable high-grade stainless steel offer unmatched lateral molding. Quick-release buckles allow for hassle-free application.

ADDITIONAL FEATURES:

- Universal Sizing - Only 2 items to stock
- Fit of a custom in less than 10 minutes
- Adjustable thumbwheel controls offloading force
- Offload force engagement from 10° - 60°
- Channeled fabric-covered padding for breathability in hot climates
- Removable stops control knee flexion in 15° increments from 15° - 60°

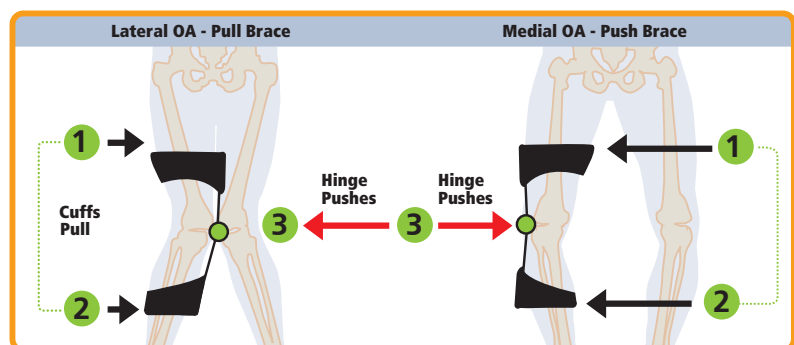


Solar OA™ Osteoarthritis Knee Brace

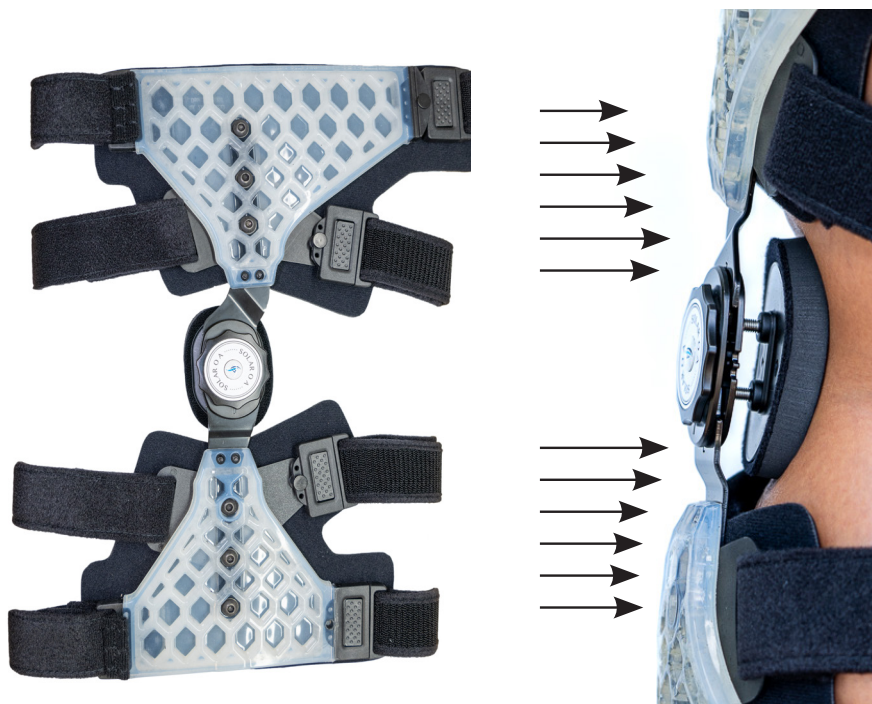

Our CMS custom-fit compliancy starts with the delivery of a brace with no initial shape that requires substantial modification to achieve an individualized fit.

Flexible silicone cuffs bonded to malleable stainless steel, offer unmatched varus/valgus molding.

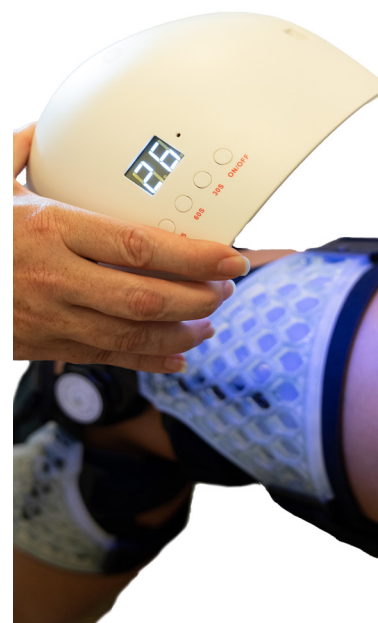

A 90-second UV light bath solidifies the brace and ensures custom definition and long-term durability.



Cotton/Lycra Undersleeve (Sold Separately)

SIZE	Product Number	Sizing Guide
XS	80-0010	13"-20"
S	80-0011	14"-22"
M	80-0012	15"-24"
L	80-0013	16"-26"
XL	80-0014	19"-28"
2XL	80-0015	22"-29"

Measurement taken 8" above mid patella

Suggested Billing Code: L-2397



Product Number

Right	80-1000
Left	80-1001

Fits maximum thigh circumference of 35", calf 25"

Suggested Billing Codes:

L-1843 (Custom Fit) or L-1851 (OTS)

Additional Components

UV Light	80-1040
Starter Kit	80-1045

Replacement Parts

Thigh Cuff - Right	80-1010
Thigh Cuff - Left	80-1011
Calf Cuff - Right	80-1012
Calf Cuff - Left	80-1013
Liner Kit - Right	80-1020
Liner Kit - Left	80-1021
Strap Kit - Universal	80-1030



CORFLEX GLOBAL | 669 EAST INDUSTRIAL PARK DRIVE | MANCHESTER, NH 03109
CORFLEXGLOBAL.COM | TOLL FREE: 800 426 7353 | FAX: 603 623 4111